



A perfect fit – ideal for those people who spend many hours a day in their wheelchair.
 Seven possible backrest positions are available for adapting to your special requirements without needing a tool. Simply release the pins that drop into the holes on the left and right side and the backrest can be adjusted as desired. The horizontal belts of the backrest covering can be adjusted to individual requirements ensuring that users have the correct physiological posture.



United Kingdom
 Gerald Simonds Healthcare Ltd
 Tel +44 (0)1296 380200
 Fax +44 (0)1296 380279
 info@gerald-simonds.co.uk

Italy
 Ortosan S.R.L.
 Fon +39 0832-230502
 Fax +39 0832-230502
 ortosansrl@tin.it

France
 Poitevin-Diffusion
 Tel +33 (0)2 47 29 12 26
 +33 (0)6 21 87 84 62
 Fax +33 (0)2 47 29 03 35
 contact@poitevin-diffusion.com

Switzerland
 PRO ACTIV Swiss GmbH
 Tel +41 (0)41 922 05 88
 Fax +41 (0)41 922 05 87
 proactiv@bluewin.ch

Denmark
 Pro Movec A/S
 Tel +45 (0)70272924
 Fax +45 (0)70272429
 info@promovec.dk

Israel
 Israel Quality of Life Center Ltd.
 Tel +972-4-90 40 780
 +972-4-90 40 880
 Fax +972-4-90 40 480
 iqlc_roy@netvision.net.il

Printed in Germany 12/2005. Technical changes and printing errors excepted.
 CE Colour deviations from actual product colours are possible due to reasons relating to printing techniques.



For further information on PRO ACTIV products and accessories, please refer to our general catalogue or contact:

PRO ACTIV Reha - Technik GmbH
 Im Hofstätt 11
 D-72359 Dotternhausen
 Telefon +49 (0)7427 9480-0
 Fax +49 (0)7427 9480-25
 info@proactiv-gmbh.de
 www.proactiv-gmbh.de

SPEEDY 4all


Rigid frame wheelchair



Optimum dimensions for stowing away – the sides disappear into the frame



Body Contour – elastic, breathable seat covering available in the same 1 cm increments as the frame width

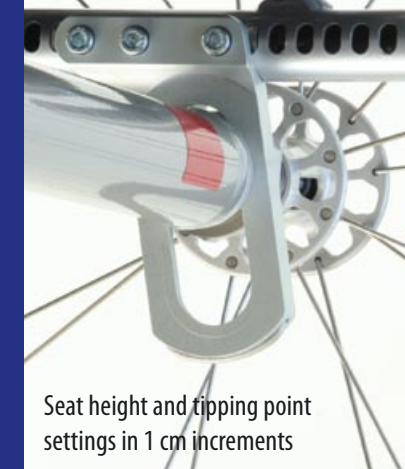

Seat height and tipping point settings in 1 cm increments



Reinforced axle mount for seat heights up to 47 cm



SPEEDY 4all – live life to the full

Unmistakable and smart, the SPEEDY 4all is a wheelchair for all lifestyles. With few components but a multitude of functions. For example, the wheel camber is integrated in the drive wheel socket, so the track can be set exactly to the millimetre simply by turning the socket. The result is a wheelchair with an overall weight (from 8.9 kg) that can compete with any specially welded frame.

Simple adjustment of the tipping point and seat height of the SPEEDY 4all make it a favourite with first timer users, whereas its light and filigree appearance convinces even the „pros“.

The continuous axle tube, high-quality craftsmanship and, last but not least, PRO ACTIV's attention to detail ensure that the SPEEDY retains its rigidity and excellent driving characteristics whatever the conditions.



Customized – e.g. special backrest of powder-coated aluminium

Suspension forks – reduce pain sensitiveness and paroxysmal extensor spasms caused by uneven terrain



Fully detachable wheel base extension with quick-release mechanism



One-arm drive

