

 **loopwheels** 

It's the only wheel you'll need

Go light
Go comfortable
Go anywhere

loopwheels.com





Lightweight, vibration-reducing wheels for everyday use.

Safe levels of whole body vibration for over 2 hours every day.

ISO SAFETY MEASURES 2631-1 1985

This measures the exposure to vibration in a seated position, before action is needed to protect health.

ATV (average daily use) tested at 12km/h over uneven surfaces

Spinergy 24" day chair wheel	safe for 1 hour 3 minutes
Softwheel 24" day chair wheel	safe for 1 hour 18 minutes (+24%)
Loopwheels Urban	safe for 2 hours 10 minutes (+106%)

Loopwheels Urban is safe for over **TWICE THE TIME** vs a Spinergy wheel on off-road surfaces at an average of 12km/h.

PRODUCT SPECIFICATION

Wheel size: 24" (540mm)

ETRTO number: 25-540 (24x1.00)

Axle bearing: 1/2" (12.7mm) and 12mm

Push rim: 16 rivnut fitting

Rivnut centres measurement: 515mm

Spring/loop materials: Carbon fibre composite

Rim: Double-walled aluminium, black

Inner rim width: 16mm

Hub: CWC machined aluminium, anodised black sealed cartridge bearings


Bearing width at faces: 1.8" (45.8mm)

Wheel weight: 1.3kg (2.87lb)

NB Wheel weights without push rim, tyre, or axle

Maximum load (user and chair): 120kg (265lb)

Recommended tyres: Schwalbe Marathon Plus

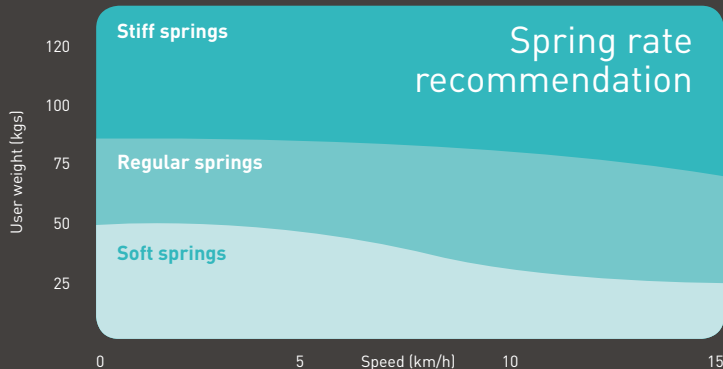

Loopwheels Urban is available in 3 spring strengths – **Soft springs**, **Regular springs** and **Stiff springs**. Use the chart above to estimate the most usual higher speed on the wheels and the fully-loaded weight of the chair.

Jelly Products Ltd

Unit 202 Boughton Industrial Estate North,
Newark, Nottinghamshire NG22 9LD

T: +44 (0) 1623 860909

E: info@loopwheels.com

W: loopwheels.com

To learn how Loopwheels can make the difference, visit loopwheels.com